

# Ozone CHEMets® Kit

**K-7404/R-7404:** 0 - 0.6 & 0.6 - 3 ppm

## Safety Information

Read SDS (available at [www.chemetrics.com](http://www.chemetrics.com)) before performing this test procedure. Wear safety glasses and protective gloves.

## Test Procedure

1. Add 5 drops of A-7400 Activator Solution to the empty sample cup (fig. 1).
2. Fill the sample cup to the 25 mL mark with the sample to be tested, being careful to minimize turbulence (fig. 2).

**Note:** Ozone loss from sample occurs rapidly. Do not transfer sample to other containers.

3. Immediately place the CHEMet ampoule, tip first, into the sample cup. Snap the tip. The ampoule will fill leaving a bubble for mixing (fig. 3).
4. To mix the ampoule, invert it several times, allowing the bubble to travel from end to end.
5. Dry the ampoule. Obtain a test result **1 minute** after snapping the tip.
6. Obtain a test result using the appropriate comparator.

- a. **Low Range Comparator (fig. 4):** Place the ampoule, flat end first, into the comparator. Hold the comparator up toward a source of light and view from the bottom. Rotate the comparator until the best color match is found.

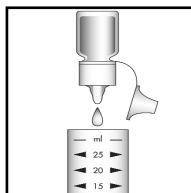


Figure 1

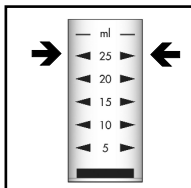


Figure 2

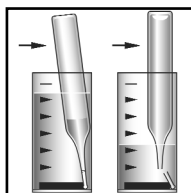


Figure 3

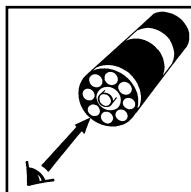


Figure 4

- b. **High Range Comparator (fig. 5):** Place the ampoule between the color standards until the best color match is found.

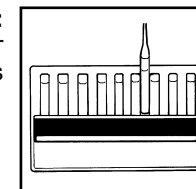


Figure 5

## Test Method

The Ozone CHEMets®<sup>1</sup> test kit employs the DPD chemistry.<sup>2,3</sup> The sample is treated with an excess of potassium iodide. Ozone oxidizes the iodide to iodine. The iodine then oxidizes DPD (N,N-diethyl-p-phenylenediamine) to form a pink colored species in direct proportion to the ozone concentration.

Various oxidizing agents such as halogens, ferric ions and cupric ions will produce high test results.

1. CHEMets is a registered trademark of CHEMetrics, Inc. U.S. Patent No. 3,634,038
2. APHA Standard Methods, 23rd ed., Method 4500-Cl G - 2000
3. EPA Methods for Chemical Analysis of Water and Wastes, Method 330.5 (1983)

Visit [www.chemetrics.com](http://www.chemetrics.com) to view product demonstration videos. Always follow the test procedure above to perform a test.



[www.chemetrics.com](http://www.chemetrics.com)  
4295 Catlett Road, Midland, VA 22728 U.S.A.  
Phone: (800) 356-3072; Fax: (540) 788-4856  
E-Mail: [orders@chemetrics.com](mailto:orders@chemetrics.com)

May 19, Rev. 4